




**KINGMED**  
»KING«-klein» address →

# *nasal mask*

Nasal Active

Nasal Fix

# nasal mask

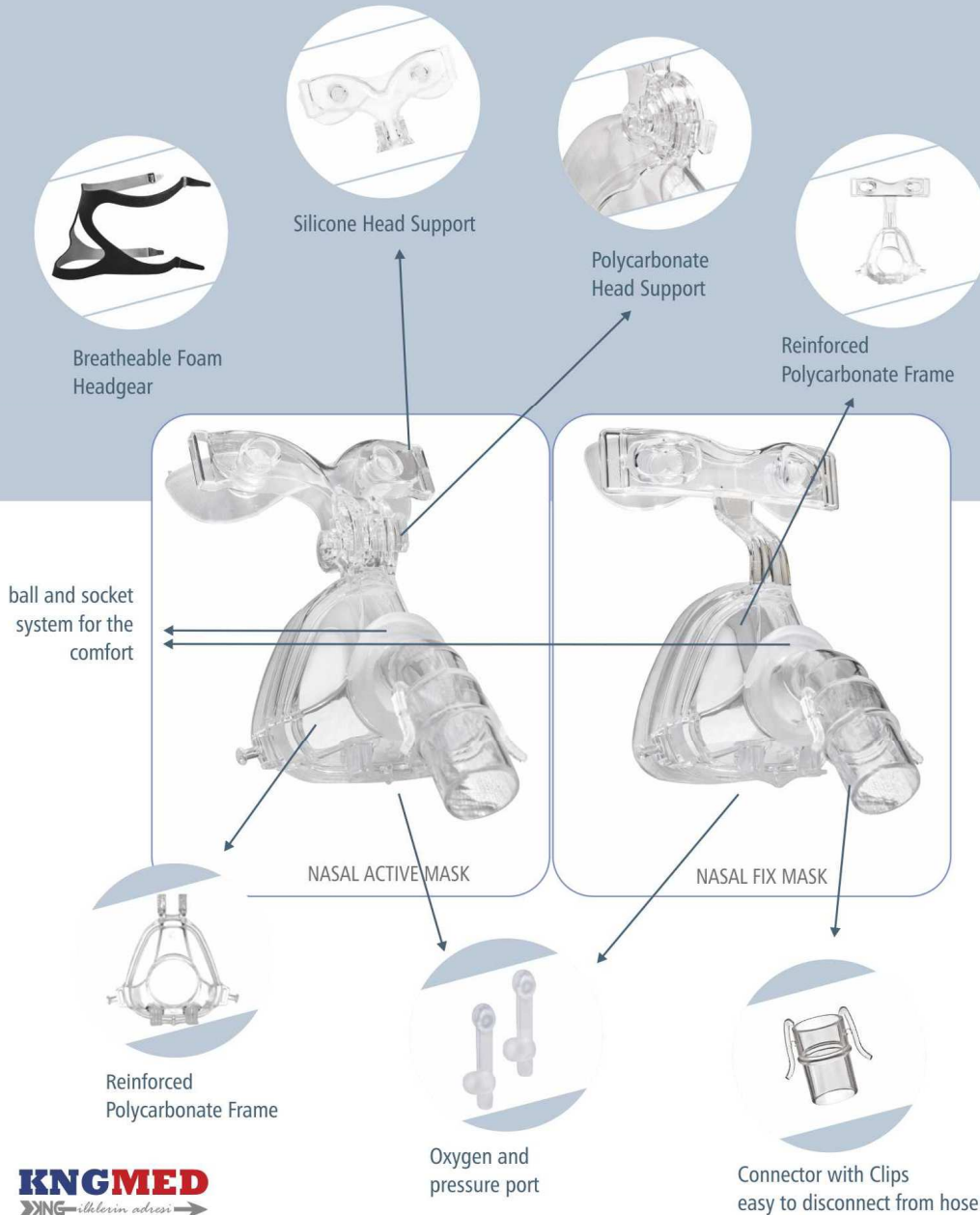

# nasal mask




## Raw Material

## Features

## Design

## Usage

## Available Sizes

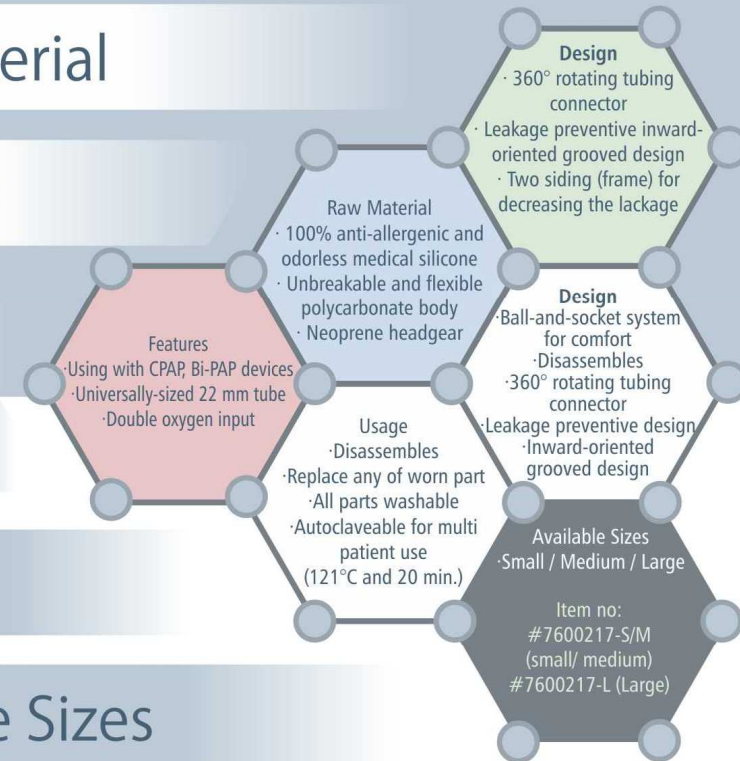

### Available Sizes (Small, Medium, Large)

- item no #7600217-S/M (small/ medium)
- item no #7600217-L (Large)





Nasal silicone mask is designed for use the patients for whom prescribed the CPAP, Bi-PAP, Auto-PAP or VPAP therapy by a health care professional or respiratory therapist. Nasal silicone mask is used with the CPAP, Bi-PAP, Auto-PAP or VPAP therapy devices for single patient use in the home or multi-patient use in the hospital or a health center.

# active - fix mask

## Cleaning the Mask at Home:

- Clean the Silicone Ora/Nasal masks and head gears before first use and daily with a mild dishwashing detergent.
- The headgear does not need to be removed for daily cleaning of the nasal masks.
- Hand wash in warm water with a mild dishwashing detergent.
- All components should be rinsed well with drinking quality water and allowed to air dry out of direct sunlight.
- Replace any parts of your mask that have visibly deteriorated as cracking, crazing, tears, etc.



## Warnings

Shelf life: 5 (five) years.

Storage temperature: - 20 °C +60 °C



- At low PAP pressures the flow through the exhalation port may not be poured out completely and rebreathing may occur.

## nasal active - fix



## For Multi-patient use in the Hospital:

- Use mask autoclaving (121°C and 20 min.) or disinfection liquids instructions to reprocess the mask between patients.



- Discontinue using mask if you have any adverse reaction, and consult your physician or sleep therapist.
- Clean the Nasal mask before first use and daily.
- Inspect all parts for damage or wear and replace any parts that have visibly deteriorated as cracking, crazing, tears, etc.



- If oxygen is used with device, the oxygen flow must be turned off when the device is not operating. Otherwise, Oxygen accumulated in the device enclosure will create a risk of fire.



- The Nasal mask and headgear should fit comfortably and the Nasal mask should not dig into your skin. Verify that the Nasal mask and headgear are the correct size.
  - Do not use the Turquoise masks, If you are using a pill that potentially causes vomiting

*Nasal*  
*active & fix mask*

**KNGMED**  
»NG-»Merkezi adresi»

Tel : +90 232 833 42 62  
Fax: +90 232 833 42 63

info@kngmed.com  
kng.info@yahoo.com

10001 Sokak No.28 / A - 10  
Çınar Sanayi Sitesi  
Ulukent / Menemen / İzmir / TÜRKİYE

