


# eeg electrodes

\* Incorporates pre-attached, tangle-free lead wires to eliminate signal artefact (due to worn connectors) and to ensure the best possible signal

\* 2 mm hole in the centre of the cup allows gel to escape during positioning and permits blunt needle insertion

\* Colour-coded lead wire allows easy tracing

\* Can be used with all gels, pastes and adhesives

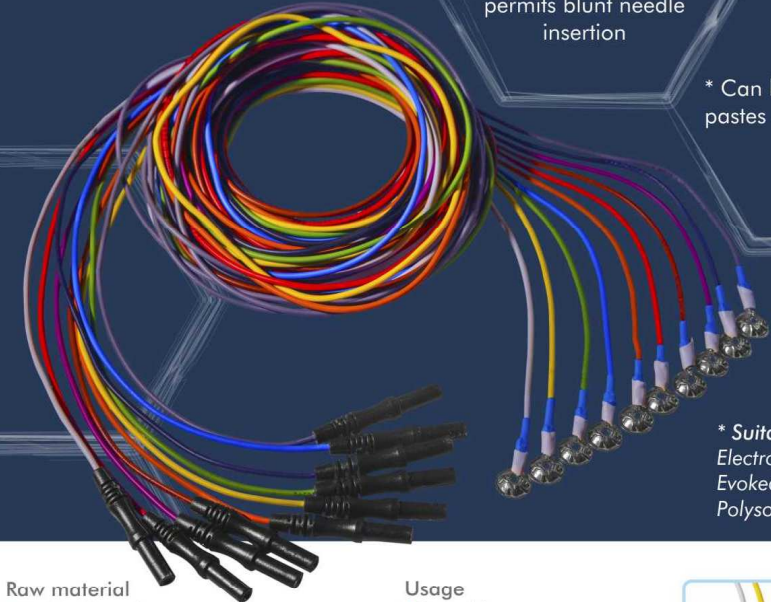

\* Suitable applications  
Electroencephalography (EE)  
Evoked Potentials (EP)  
Polysomnography (PSG)

#### Raw material

- ▶ Gold, silver disc
- ▶ Strengthened cable with silk

#### Design

- ▶ Shrink tube technology
- ▶ Tangle free form
- ▶ Female/male connector
- ▶ Designed to eliminate excessive paste

#### Features

- ▶ Insulated coating to prevent interference
- ▶ 9-10mm disk diameter
- ▶ 10 different colors

#### Usage

- ▶ Reusable

#### Size

- ▶ Custom length

#### Item no

- ▶ #7600019-Au (Gold)
- ▶ #7600019-Ag (Silver)



eeg electrodes

